

Setting Up Your Glucometer

Going through glucometer menu:

Press and hold OK button > Summary > Time/Date

Results log

Language

Trend 90

Range/90d >

before meal (3.8-5.2 mmol/L)

Settings

Features

after meal (4.0-7.7 mmol/L)

Messages

(don't worry about the 90d target)

Bluetooth

Control test

Meter Information

Time/Date

Language

Range/90d

Features

> Basic Mode **OFF**

Messages

Test Tracker **4/day**

Bluetooth

Range indicator **Emoji**

Control test

Meter information

Time/Date

Language

Range/90d

Features

Messages

> Pattern **OFF**

Bluetooth

Mentor **OFF**

Control test

Awards **OFF**

Meter information

Getting Started Guide

OneTouch® Delica® Plus lancing device

Make every test Virtually Pain Free



Preparing your lancing device



1 Remove the lancing device cap



2 Insert a sterile lancet into the lancing device and remove its protective cover



3 Replace the lancing device cap



4 Adjust the depth setting
Smaller numbers are for a shallower puncture and larger numbers for a deeper puncture.



5 Cock the lancing device by sliding the button back until it clicks

Now you can test

Testing and removing the used lancet



Puncture your finger

Hold the lancing device firmly against the side of your finger. Press the release button to puncture.



Apply the sample

Apply the blood sample following the instructions in your OneTouch® Owner's Booklet for proper sample application.



3 Remove the lancing device cap again



Cover the exposed lancet tip

Place the lancet protection cover on a hard surface and push the lancet tip into the flat side of the disk.



Eject the lancet

Holding the lancing device downwards, push the slider forward until the lancet comes out. Push the cap back onto the device.



6 Replace the lancing device cap

Discard the used lancet carefully after each use to avoid unintended lancet stick injuries.

Use with OneTouch® Delica® Plus or OneTouch® Delica® lancets.

Before using the OneTouch® Delica® Plus lancing device to test blood glucose, carefully read the OneTouch® meter Owner's Booklet and inserts that come with the system.

ONETOUCH®
every touch is a step forward®


Tips to help make blood glucose checks less painful and reduce the chance of infection

- **Before testing, wash your hands and the sample site with warm, soapy water.** Rinse and dry completely. Contaminants on the skin may affect results.
- **Use a new sterile lancet each time you test. Do Not Re-use**
- **Gently squeeze and/or massage your fingertip until a round drop of blood forms.**
- **Apply the lancing device to the side of the upper part of your finger.**
- **Choose a different puncture site each time you test.** Repeated punctures in the same spot may cause soreness and calluses.
- **Never share a lancet or lancing device with anyone.** For single patient use only.




Use a new OneTouch® Delica® Plus sterile lancet with each test to make lancing virtually pain free.

Lancet reuse can dull or bend the tip of the lancet, causing damage to your skin, scarring and greater pain.

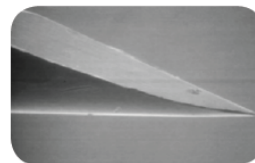

Native lancet before use

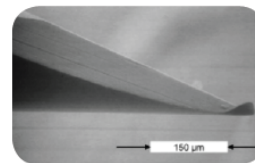

Used lancet after a single skin penetration¹

These images are from a publication that provides an overview of lancing and are not representative of any OneTouch® device.

Ask for the OneTouch® Delica® Plus lancing device and 30g (Fine) or 33g (Extra Fine) lancets at your pharmacy.



1. Heinemann L, Boecker D. Lancing: Quo Vadis? Journal of Diabetes Science and Technology 2011;5(4):966–981.
© 2019 LifeScan IP Holdings, LLC US-DPL-1900002

ONETOUCH®
every touch is a step forward[®]

Getting to know your OneTouch Verio Reflect® meter

Before using the OneTouch Verio Reflect® meter to test your blood glucose, carefully read the Owner's Booklet, and the inserts that come with the components of the system.



- 1 Test strip port
- 2 Bluetooth icon
- 3 Unit of measure (mmol/L)
- 4 Micro USB data port
- 5 OK button
- 6 Back button
- 7 Up and down buttons
- 8 ColourSure® Dynamic Range Indicator
- 9 Range Indicator icon
- 10 Test Tracker icons

Quick start: Setting up your meter

1 Turn your meter on by pressing and holding OK until the **Welcome** screen appears. You can also turn your meter on by inserting a test strip.

2 OPTIONAL: You may choose to set up your meter by using the OneTouch Reveal® App.

3 Set your language, time and date by pressing the **Up and Down Buttons** to highlight, and OK to save.

4 OPTIONAL: You may also set your Range Limits, Daily Test Goal and 90 Day Average Target* at this time.

5 **Setup complete** appears. Your meter is now ready for use.

6 Turning your meter off can be done any 1 of 3 ways: by pressing and holding OK for several seconds, removing a test strip, or leaving the meter alone for 2 minutes.

* Be sure to talk to your healthcare professional about the settings that are right for you.

Quick start: Taking a test

1 Wash your hands with warm soapy water and dry them thoroughly.

2 Insert a test strip into the **test strip port**, ensuring the gold side and two silver prongs are facing you.

3 Adjust your lancing device's depth setting, then hold the device firmly against the side of your finger.*

* Always use a new, sterile lancet each time you test.

4 Press **release button** to puncture your finger. Remove lancing device. Gently squeeze your finger to form a round droplet of blood.

5 Hold the meter at a slight angle, with either side of the test strip lightly touching the blood droplet. The test strip draws blood into the channel.

6 Your meter instantly lets you know if your current result is below, within or above your range limits. It also shows when your result is near the low or high end of your target range.

Visit OneTouch.ca or contact OneTouch® Customer Care


- Visit OneTouch.ca
- Email contact@onetouch.ca
- Phone 1 800 663-5521

Apple, App Store, and the Apple logo are trademarks of Apple Inc., registered in the US and other countries. Google Play and the Google Play logo are trademarks of Google LLC. The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by LifeScan Scotland Ltd. and its affiliates is under license. Other trademarks and trade names are those of their respective owners. Always read and follow the label. ©2019-2022 LifeScan IP Holdings, LLC. All rights reserved. CA-VRF-1900170

ONETOUCH

ColourSure® PLUS

ColourSure® Dynamic Range Indicator feature instantly lets you know where your result falls within your range limits.

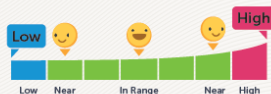

Blood Sugar Mentor™

Blood Sugar Mentor™ messages automatically appear after a test result is displayed.



ONETOUCH Reveal®

The OneTouch Reveal® mobile app gives you even more insight to help you manage your diabetes.*



Pairing instructions

The Bluetooth® feature must be turned on in order to pair the meter with the app.

Open the OneTouch Reveal® mobile app and follow instructions to pair the meter with your compatible wireless device.

* More than just using the meter

Connect for even more insight



Download OneTouch Reveal® for free:



Then visit OneTouch.ca today

* Data on File.

